

the difference is clear

# Vitalograph AIM

## For effective inhaler training on dry powder and metered dose inhalers

The unique Vitalograph AIM™ (Aerosol Inhalation Monitor) is designed to assist medical practitioners in training subjects to use their inhaler correctly. Mastery of the correct technique will support more accurate drug delivery and good patient compliance, resulting in better disease management and fewer visits to medical professionals.

### AIM validates good technique in the use of DPI and MDI inhalers:

The AIM uses a hygienic single-use disposable inhaler simulator to provide clear feedback and information to support correct technique regarding:

- Inspiratory acceleration at start of inspiration
- Timing of firing of MDI / activation of DPI
- Inspiratory flow rate throughout inspiration
- Inhalation time within target flow range
- Breath hold time at end of inhalation

### Benefits of the AIM:

- Validates effective use of an inhaler
- Identifies correct choice of inhaler device type
- Identifies and quantifies poor inhaler technique
- Objective identification of proper use of the inhaler device
- Assists in training patients to use their inhalers correctly
- Provides clear feedback and instructions to subject and trainer
- Easy to use, hygienic and low cost
- Can optionally be used with a spacer

**Vitalograph is a world leading provider of outstanding quality cardio-respiratory diagnostic devices, clinical trial services and medical equipment servicing. With a pioneering heritage of excellence spanning half a century Vitalograph continues to make valuable contributions to effective medical care and enhanced quality of life.**




Simple to use, well designed and attractively priced, the AIM ensures that the correct inhaler device has been selected and is being used as effectively as possible.

A summary result is displayed to clearly illustrate the difference between good and poor technique. The parts of the inhalation process which a subject has yet to master, including factors that are difficult to observe by eye, are clearly highlighted thereby enabling physicians to concentrate on specific aspects of subject technique.

- Displays aspects of the technique that were good or poor

### Technique Good/Poor Summary:

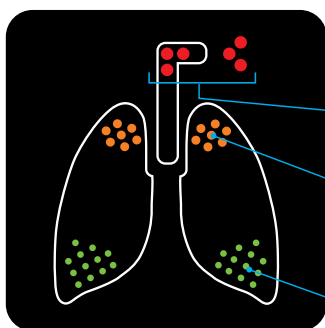

#### Fail (Red):

Sub-optimal inspiratory flow rate or the inhaler was not activated/fired

#### Poor (Orange):

Not met inhalation time or breath hold too short

#### Good (Green):

Correct force of inhalation with adequate inspired volume and breath hold time

## Technical Specification:

**Product:** AIM (Aerosol Inhalation Monitor)

**Model Number:** 4500

**Power Supply:** Battery AAA x 4

**Weight:** Net: 260g (inc batteries and tubing)  
Gross: 310g

#### Dimensions:

Net: 165mm x 133mm x 39.6mm

Gross: 190mm x 140mm x 65mm

**Operating Temperature:** 10°C – 40°C

**Storage Temperature:** 0°C – 50°C

**Storage Relative Humidity:** 10% - 95%

**Safety Standards:** EN 60601-1:2006

**Medical Safety Standard:** Medical Devices Directive 93/42/EEC1993 as amended

#### Designed & Manufactured Under:

ISO 13485:2003, FDA 21CFR820, CMDR & JPAL

## Ordering Info:

45710 4500 AIM (Aerosol Inhalation Monitor)

45610 Disposable DPI Inhaler Simulator (25)

45611 Disposable MDI Inhaler Simulator (25)

79192 Replacement silicone tubing

45027 HFA Placebo Aerosol (8)

The AIM uses hygienic single-use disposable inhaler simulators



Validates effective inhaler use



Your respiratory partner

#### UK Sales

**Vitalograph Ltd.**

Maids Moreton Buckingham

MK18 1SW England

Tel: 01280 827110

Fax: 01280 823302

e-mail: [sales@vitalograph.co.uk](mailto:sales@vitalograph.co.uk)

[www.vitalograph.co.uk](http://www.vitalograph.co.uk)

#### Export Sales

**Vitalograph Ltd.**

Maids Moreton Buckingham

MK18 1SW England

Tel: +44 1280 827120

Fax: +44 1280 823302

e-mail: [sales@vitalograph.co.uk](mailto:sales@vitalograph.co.uk)

[www.vitalograph.eu](http://www.vitalograph.eu)

#### Vitalograph GmbH

Rellinger Straße 64a

20257 Hamburg

Germany

Tel: (040) 547391-0

Fax: (040) 547391-40

e-mail: [info@vitalograph.de](mailto:info@vitalograph.de)

[www.vitalograph.de](http://www.vitalograph.de)

#### Vitalograph Inc.

13310 West 99th Street

Lenexa

Kansas 66215 USA

Tel: (913) 888 4221

Fax: (913) 888 4259

e-mail: [vitcs@vitalograph.com](mailto:vitcs@vitalograph.com)

[www.vitalograph.com](http://www.vitalograph.com)

#### Vitalograph (Ireland) Ltd.

Gort Road Business Park

Ennis

Co Clare Ireland

Tel: (065) 6864100

Fax: (065) 6829289

e-mail: [sales@vitalograph.ie](mailto:sales@vitalograph.ie)

[www.vitalograph.ie](http://www.vitalograph.ie)

**[www.vitalograph.co.uk](http://www.vitalograph.co.uk)**

Vitalograph® and AIM™ are trademarks or registered trademarks of Vitalograph Ltd.



0086

PRINT REF: 14697\_4